

MultiShape

Powerful Cavitation, Multi-Polar RF and Vacuum

The most optimal solution for body
and skin care Innovative 3-in-1 System

Dynamic vibration

Dramatic Eye Lifting

Intensive Lipoplasty

SuperPulse Technology




What is MultiShape?

- The best combinations of Multi-polar RF, Vacuum and Cavitation have shown effective results in cellulite reduction and body reshaping
- Using the specially designed Multi-polar applicator, target management for eye & neck lifting
- Soft and powerful vibrational effect through the SuperPulse suction system
- Multipurpose management to accommodate for the diverse needs of patients
- Using a diverse array of applicators, it is an easy and convenient way for whole body management (removable method and automatic recognition function)
- Technology integrated to reduce vibrational sounds to bare minimum.

Technology


Superpulse

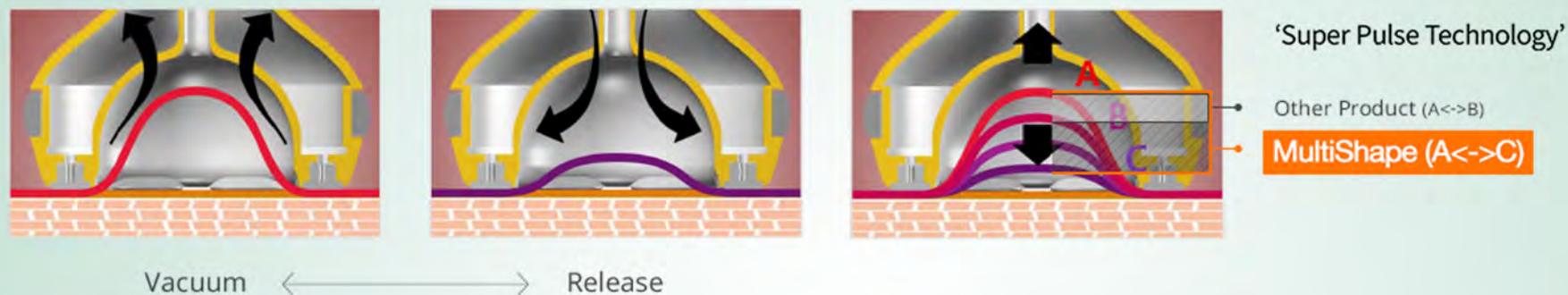

- Using the new patented technology SuperPulse from the Vacuum, the best results are achieved for cellulite reduction and body reshaping management
- Endermologie function with no pain and high quantity of motion
- Able to control the level of Vacuum and Release to achieve the best results
- User is able to enjoy the smooth and gentle sensation of the *Sliding Technology
- 15 pulse outputs per second to achieve maximum effect of physical activity

*Sliding Technology - Enables the selection of vacuum and release area between "A" , "B" , "C" , facilitating the smooth and free movement of the handpiece which enhances the treatment efficiency

Application

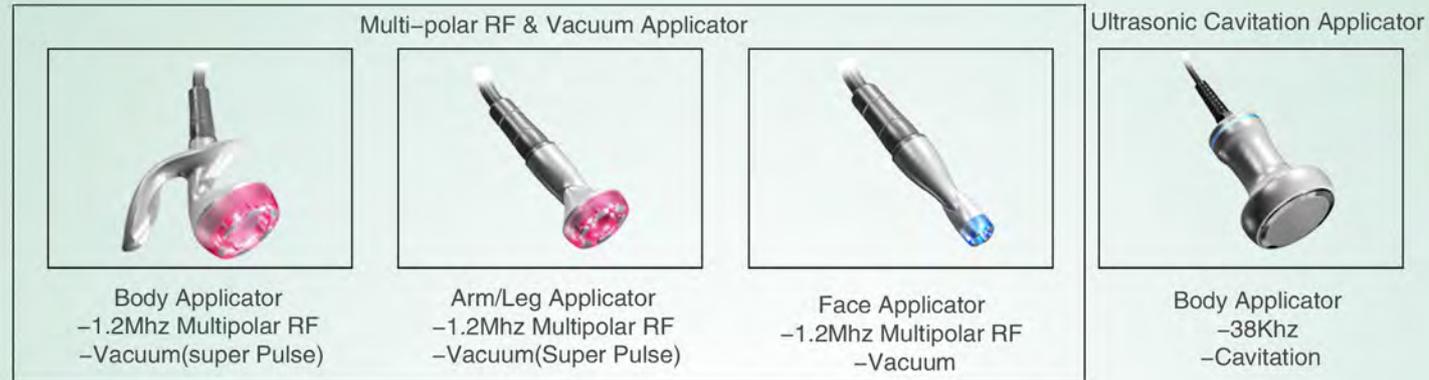

Effects of Multi-polar RF

- Increases intracellular diffusion of oxygen by heat
- Increases the microcirculation
- Actively promotes the degradation and lipid metabolism of fat cells
- Ensures tightening of the collagen fibers
- Improves the elasticity and firming of the skin
- Light Therapy effect via RF technology
- No pain, effective & completely safe

Effects of Vacuum

- Promotes lymphatic and degraded fat cells drainage and blood circulation
- Reduces the subcutaneous layer of the fat chamber
- Increases blood flow in the veins and facilitate the supply of oxygen to the blood
- Reduces the adhesiveness of fat cells
- Relaxes the tissue and massages the fascia
- Relieves swelling and greatly improves body and facial contour line
- No pain, effective & completely safe

Effects of Ultrasonic Cavitation

- 38KHz ultrasonic cavitation' s wavelength promotes degradation of fat by reaching deep into the area of subcutaneous fat normally difficult to be reached by average high frequency procedure
- No pain, effective & completely safe

Other features



Large(Body) applicator : red LED 630 nm
Middle (Body & Face) applicator : red LED 630 nm
Small (Face) applicator : blue LED 430 nm
Need a safety key to operate the system



Need a safety key to operate the system



Sliding Touch Power ON / OFF
Intuitive Touch Control
Convenient Auto-Sensing Connection



Start or stop the system sith the switch (convenient use)

Structure



LCD Touch Screen

Safty Key



Body Applicator(Ultrasonic Cavitation)



Face Applicator(RF/Vacuum)



Arm/Leg Applicator(RF/SuperPulse/Vacuum)



Body Applicator(RF/SuperPulse/Vacuum)

Specifications

Power Input	220~240VAC 50/60 Hz
Maximum Power Consumption	900W
Classification (EN60601-1)	Class1 , Type BF
User Inter Face	10.2Inch Touch LCD
Size (W D H)	388 x 354 x 1095 mm
Weight	50 Kg
RF	Frequency : 1.2MHz Bi-Polar , Multi Polar , Max 50W Light Therapy : Red LED (630nm) , Blue LED (430nm)
Cavitation	Frequency : 36KHz \pm 10% Max Power : 50W Pulse Type : 4 (Continues , Pulse1 , Pulse2 , Auto)
Vacuum	Max Vacuum : 350mmHg 1 ~ 15Hz Pulse (Continues , Pulse , Super-Pulse)
Treatment Timer	0~30min (Preset 10min)